


# Catalyst 5 Cross Braces - Open Seating (Seating System)

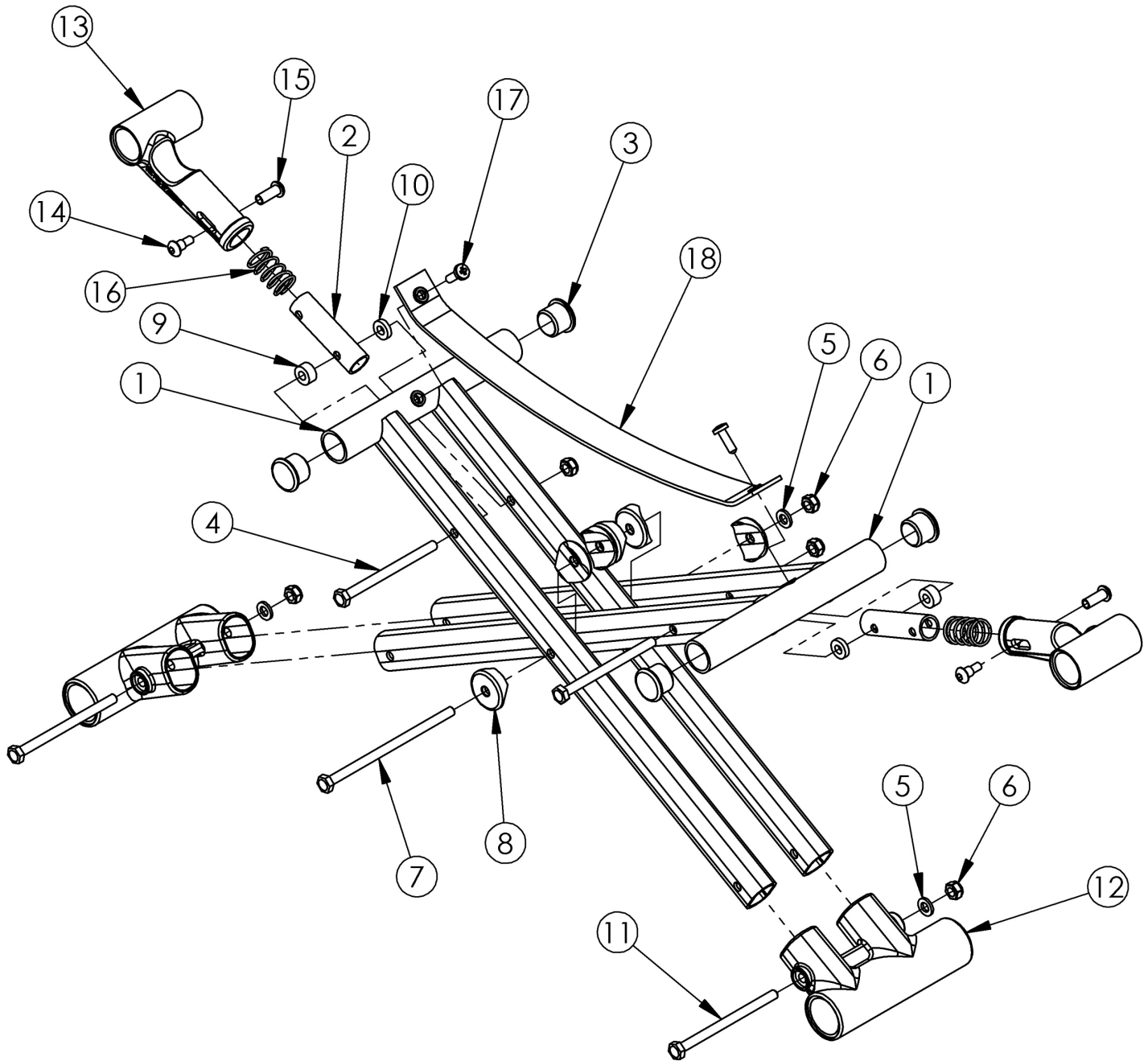

*Pictured oval-style cross braces effective 6/15/16 on 14" - 22" wide chairs starting with serial numbers C5032544, C5037980, C5HD002834 and CVX0091137. When changing from round-style to oval-style cross braces, lower frame pivots and cross brace hardware are required. Order kit 502177.*

Item Number	Part Number	Description	Quantity
2a, 13a, 14-16	500024	Upper Pivot Assembly 14" Width	2
2b, 13b, 14-16	500025	Upper Pivot Assembly 15" Width	2

Item Number	Part Number	Description	Quantity
2c, 13b, 14-16	500026	Upper Pivot Assembly 16" Width	2
2d, 13b, 14-16	500027	Upper Pivot Assembly 17" Width	2
2e, 13b, 14-16	500028	Upper Pivot Assembly 18" Width	2
2f, 13b, 14-16	500029	Upper Pivot Assembly 19" Width	2
2g, 13b, 14-16	500030	Upper Pivot Assembly 20" Width	2
2h, 13b, 14-16	501479	Upper Pivot Assembly 21/22" Width	2
4-8, 11, 12	107296	Lower Frame Pivot and Catalyst Oval-Style Cross Brace Hardware Kit 14" Wide (Pair)	1
4-12	502177	Lower Frame Pivot and Catalyst Oval-Style Cross Brace Hardware Kit (Pair)	1
1a, 3	003353	Seating System Cross Brace 14" Wide CAT	1
1b, 3	003354	Seating System Cross Brace 15" Wide CAT	1
1c, 3	003355	Seating System Cross Brace 16" Wide CAT	1
1d, 3	003356	Seating System Cross Brace 17" Wide CAT	1
1f, 3	003357	Seating System Cross Brace 18" Wide CAT	1
1g, 3	003358	Seating System Cross Brace 19" Wide CAT	1
1h, 3	003359	Seating System Cross Brace 20" Wide CAT	1
1i, 3	003360	Seating System Cross Brace 21" Wide CAT	1
1j, 3	003361	Seating System Cross Brace 22" Wide CAT	1
3	101162	7/8 OD End Plug BLK	2
4	100660	M6 x 80mm CLS8.8 DIN931 HHCS ZC	2
5	100746	M6 Flat Washer Blk Zn	3
6	100658	M6 Nylock Nut Blk Zn	5
7	100678	M6 x 110 CLS8.8 DIN931 HHCS ZC	1
8	002807	Cross Tube Saddle CAT	6
9	100880	1/4" Plastic Spacer <i>Not used on 14" wide chairs</i>	2
10	105019	3/16" Plastic Spacer <i>Not used on 14" wide chairs</i>	2
11	100704	M6 x 90 CLS8.8 DIN931 HHCS ZC	2

Item Number	Part Number	Description	Quantity
4-11	503320	Catalyst Cross Brace Hardware Kit	1
12	002806	Lower Frame Pivot CAT	1
17	100734	M5 x 16 PHMS ZC	2
18a	200233	Seat Strap Open Frame 14" Wide CAT	1
18b	200234	Seat Strap Open Frame 15" Wide CAT	1
18c	200235	Seat Strap Open Frame 16/21" Wide CAT	1
18d	200236	Seat Strap Open Frame 17/22" Wide CAT	1
18e	200237	Seat Strap Open Frame 18" Wide CAT	1
18f	200238	Seat Strap Open Frame 19" Wide CAT	1
18g	200239	Seat Strap Open Frame 20" Wide CAT	1