


Discover What's NEW at Sunrise Medical

2 0 1 4



New!

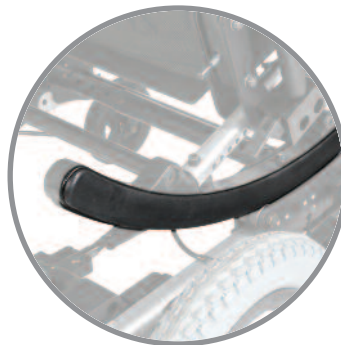

improving people's lives

Quickie® 2 and Quickie® 2 Lite Enhancements



SWING-IN/SWING-OUT HANGER

- Optimized bi-directional handle design and placement for ease of release
- Tapered locking mechanism ensures a secure lock for the life of the product



Watch VIDEO to see new features in action!




ELR/ALR

- 2" of height adjustment to better align pivot with knee joint
- New angle adjustable mechanism is easy to release and indexed with audible clicks through range of motion



24" 5-SPOKE MAG WHEEL

- Over 1 pound lighter than previous model (pair)
- Closed spoke design prevents debris from collecting on spoke

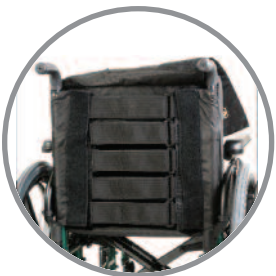

TENSION-ADJUSTABLE UPHOLSTERY OPTIONS

- Now HCPCS Code Verified
- JAY GO Back (Shown), Ballistic and 3DX Options



UPDATED HEMI FRAME HANGER

- STFHS are approximately 1/2" lower
- Now in line with seat rail



Quickie 2



Quickie 2 Lite

WC-19 TRANSIT APPROVED FOR OCCUPANT TRANSPORT



Zippie® and Quickie® IRIS™ Enhancements



SWING-IN/SWING-OUT HANGER

- Optimized bi-directional lever design and placement for ease of release
- Tapered locking mechanism ensures a secure lock for the life of the product



Watch VIDEO to see new features in action!




ELR/ALR

- 2" of height adjustment to better align pivot with knee joint
- New angle adjustable mechanism is easy to release and indexed with audible clicks through range of motion



STRUT TUBE TECHNOLOGY

- Standard availability of light weight fixed width option
- Optional width-adjustable system accommodates changes in width based on client needs

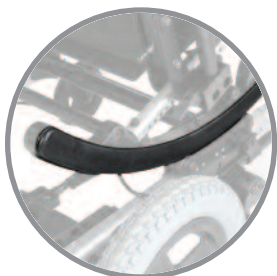

ROCKER ARM TECHNOLOGY

- Cover to prevent debris from gathering on arm
- Enhanced design to minimize lateral movement



**WC-19 TRANSIT
APPROVED
FOR OCCUPANT
TRANSPORT**

2014
Cover

Quickie® QM-7 Series Enhancements



ASL® Mushroom Joystick



ASL® Compact Joystick with Integrated Buttons




Switch-It® 3 switch Head Array



ASL® Proportional Micro Extremity Joystick



Switch-It® Buddy Button

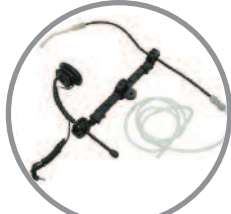

Quickie Sip N' Puff System



Switch-It® Proximity Switches



ASL® Egg Switch



Active Controls Power Buddy



RAM® Phone and Tablet Mounts



Active Controls ReCharge™ with Socket



American Flag Shroud Color

R-net PM120 Controller

Standard PM120 Controller for non-expandable or expandable controls offers 30% more power and an improved drive algorithm.

Specialty Input Devices

- Industry proven ASL® & Switch-It® technologies can be ordered on chair and will ship complete with base and ready for user set-up
- View a Video Series of our Proportional and Non-Proportional Controls by scanning this code



Watch VIDEO to see new features in action!

Proportional controls (also known as analog) allow users to move in any direction. The movement of the chair is literally proportional to the movement of the input device. Utilization of a proportional control requires coordination and precise motor control. All proportional alternative input devices require the use of an OMNI.

- 6 new Switch-it® options
- 6 new ASL® options
- New Quickie® Proportional Chin Joystick

Non-proportional (also known as digital) controls require less coordination and motor control. Each direction and speed is pre-programmed. All non-proportional alternative input devices require the use of an Omni.

- 9 new Switch-it® options now available

Consumer Convenience Accessories

QM-7 Series products now offer a series of user accessories that help users adapt more easily to today's consumer electronic world.

- RAM® Phone and Tablet Mounts: Connect and disconnect mobile phones or tables with ease
- Active Controls Power Buddy: Designed to provide an optional LED light as well as an auxiliary USB power source
- RAM® self-leveling cup holder: Gimbal system keeps drinks level while traveling or parking on inclines
- Active Controls ReCharge™ with Socket: Helps users connect and disconnect their battery charger with ease

New American Flag Shroud Color