



RAZ REHAB SHOWER COMMUNE CHAIRS

Designed and manufactured to the highest standards of the discerning healthcare professional.



RAZ-AP



RAZ-SP



RAZ-AT



RAZ-AP600



RAZ-SP600

NEW HEAVY-DUTY MODELS!
600 lb. Weight Capacity

RAZ has created the most configurable, adjustable and comfortable rehab shower commode chairs on the market today. The chairs were developed with meticulous attention to detail – from the electro-polished stainless steel frames with easy access for users and caregivers, to precisely contoured seats molded with soft integral skin foam. Only Raz offers a commode seat that adjusts to accommodate pelvic rotation with a centering of the aperture around the user's ischial tuberosities.

Raz Rehab Shower Commode Chairs, which provide positioning, pressure reduction and promote function, are ideal for the rehab client and the institutional user.

Models shown with optional accessories: RAZ-AP (Attendant Propelled), RAZ-SP (Self-Propelled), RAZ-AT (Attendant Tilt), RAZ-AP600 & RAZ-SP600 (Heavy-Duty)


RAZ REHAB SHOWER COMMODE CHAIR COMMON FEATURES

Ischial Pelvic Alignment Seat



- Contoured for pressure reduction
- Soft integral-skin foam for comfort
- Designed to help stabilize the pelvis
- Available with open or bridged front

Adjustable / Quick-Release Seat

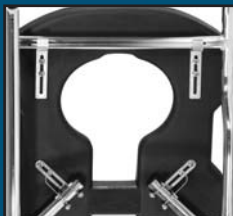


- Provides 2" of seat depth adjustment
- Locks in place for user transfers
- Removes easily for cleaning
- Optional spacers for seat angle adjustment

Adjustable Tension Backrest



- Provides trunk support for positioning/comfort
- Removable for machine washing
- Head support option adjusts for height, fore/aft position and headrest angle

Telescoping Caster Extensions

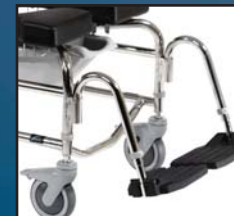

- Provide 4" of seat height adjustment without the use of tools
- Dropped front frame tube improves access over toilet bowls
- Optional removable swing-away footrests

PRx - Pressure Reduction Seat



- Dual foam construction
- Extra soft integral-skin "gel" foam
- Contoured foam insert
- 10 configurations avail.: bridge front, right/left/ front /rear openings (all in 16" and 19" depths)

Elongated Commode Pan



- Handle and spout for easy emptying
- High-gloss / easy-to-clean surface
- Mounts directly under seat frame for better containment
- Pan secures in place

Swing-Away Lateral Support



- Provides positive thoracic support
- Multi-axis adjustment
- Swings completely away for transfers
- 6" x 5" molded pad
- Available extension plate (incr. offset 1½")

Directional Caster Lock



- Keeps rear caster in a forward position
- Improves straight-line steering (e.g. down hallways or corridors)
- Allows dual-lock caster to lock independently for maximum stability

EXCLUSIVE RAZ AP/SP FEATURES




- Adjusta-Back (Option)**
- Adjusts to 85°, 90°, and 95° back angle to suit user's needs
 - Armrests adjust in height / angle
 - Includes locking armrest feature
 - Available in 17" and 20" back height

EXCLUSIVE RAZ AT FEATURES




- Adjustable Pivot Posts**
- Provide 4" of seat height adjustment without frame-toilet bowl interference
 - Pressurized gas springs work against gravity making tilting up to 40° easier for an attendant

EXCLUSIVE RAZ AP600/SP600 FEATURES

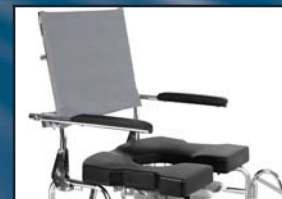


- Multi-Position Back**
- Adjusts to 95° and 100° back angle for personalized comfort
 - Flip-up/ fixed height armrests
 - Armrests adjust in angle based on user's unique needs
 - Available with 20" back height

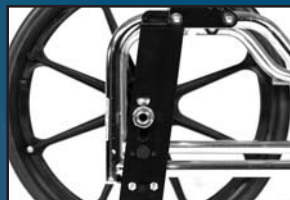

- Multi-Position Axle Plates**
- Provides up to 4" of vertical seat height and 3½" of horizontal adjustment for optimal user-toilet bowl access
 - Allows for the use of standard 22" or optional 24" rear wheels



- Tilt Assist Pedal (Option)**
- Enables attendant to use additional strength from legs to tilt larger/heavier users
 - Solid stainless steel construction fits directly into chair frame
 - Anti-slip tread on pedal

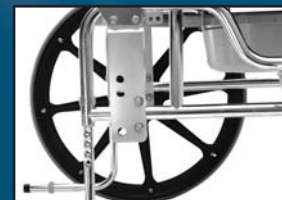

- C.O.G. Adjustable Frame**
- Provides independent horizontal adjustment of back (1") and axle plate (1½") for optimal maneuverability and stability
 - Axle plates provide up to 3¾" of vertical seat height adjustment

MANUFACTURED BY:

Raz Design Inc.

19 Railside Rd., Toronto, Ontario

Canada, M3A 1B2

Phone: (877) 720-5678

Fax: (877) 720-6586

Email: info@razdesigninc.com

Website: www.razdesigninc.com

AVAILABLE FROM:

Two year limited warranty on the frame against manufacturing defects. One year on all other components subject to normal wear and tear exclusions.

For specifications and other inquiries please visit www.razdesigninc.com