



# StarLock

The Vertical Air Cell Cushion that locks in place



## Vertical Air Cells

Star cushions feature vertical air cell technology designed specifically for individuals at high risk of pressure injuries.

Star's interconnected vertical air cell design ensures internal cell pressures equalise as the user sits in the cushion. The vertical air cell construction also effectively reduces surface tension and enables excellent envelopment.



## Precision Manufacturing

The Star range is manufactured in Sweden with cutting-edge robot production to ensure precision and consistency.

# StarLock – No Compromises

In air cell cushions, the combination of pressure redistributing qualities and stability have always been a difficult compromise – With StarLock, you can have both.

## Pressure Care with Pressure Redistribution and Postural Support

StarLock is the complete cushion solution for pressure care with pressure redistribution, stability and postural support. It is available in 5, 7.5, 10 and 13 cm cell height.

What ever the cell height, the sitting height will always be the same. However, the higher the cell, the greater the pressure redistribution.

The 13 cm cell height is unique to StarLock and offers greater immersion and envelopment.

## Shaping, Offloading and Stability


The StarLock, with its StarLock Technology, provides the possibility to lock air into each vertical air cell which makes it possible to create a desired shape, offload localised pressure and increase stability, allowing for individualised postural support.

It can be adjusted and reset again and again, making it very flexible over time.

StarLock can also work as a single air chamber cushion allowing individual setup for optimal pressure redistribution.




## Features

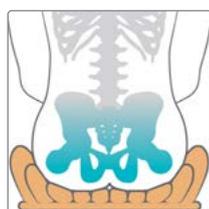


### Pressure Redistribution

StarLock is the only air cell cushion with a 13 cm cell height - that redistributes pressure over a larger surface, provides greater immersion into the cushion, and continually readjusts to the user's shape as they move in their seat.

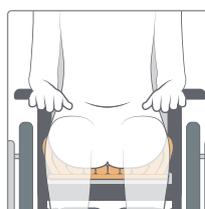


### Adjustability

StarLock can be adjusted and reset again and again, making it very flexible over time. The user can remain seated throughout the assessment.

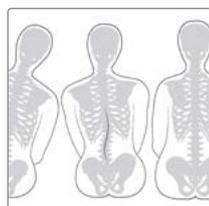


### Shape for the Individual

StarLock technology conforms to the user's body shape or desired position, which can help by locking the air into each cell. StarLock can meet the individual user's changing needs for support and shape.

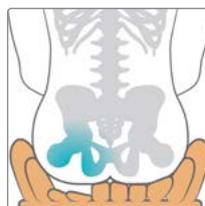

### Offloading

Locking down cells in a selected area of the cushion can relieve pressure on particularly vulnerable areas.

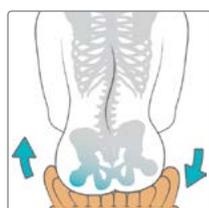

### Stability

The individual cell locking technology aids stability and allows individualised postural support.




### Recondition and Reuse

StarLock can always be washed and repaired, and it can be adjusted and disinfected as often as you want. Without locking the Lock Pockets, StarLock can also be used as a classic single chamber cushion.

STARLOCK  
TECHNOLOGY

Air Cells

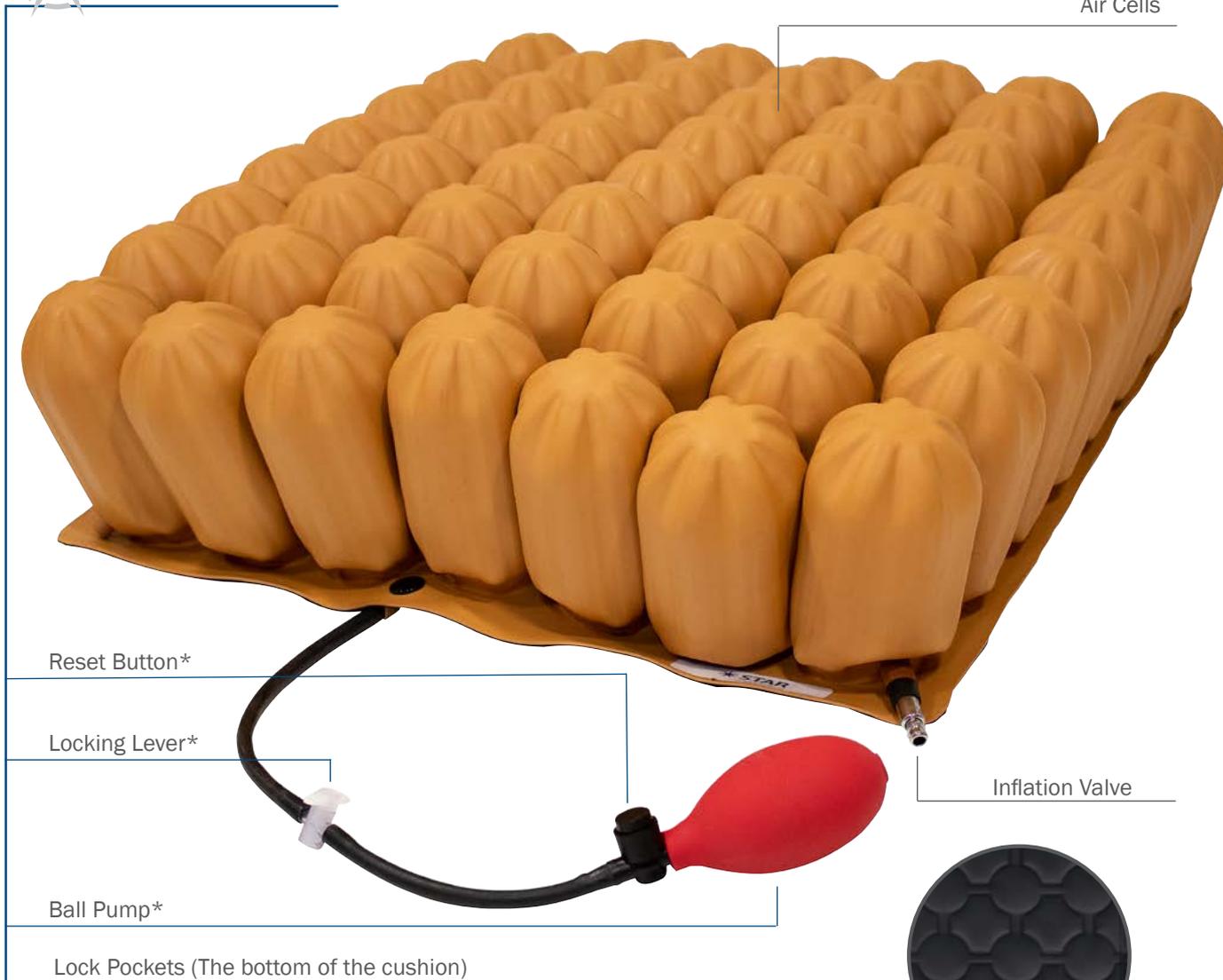

Reset Button\*

Locking Lever\*

Ball Pump\*

Lock Pockets (The bottom of the cushion)

Inflation Valve



\*Reset Button, Locking lever and Ball Pump are all for activate and deactivate the locking of the Lock Pockets.

StarLock cushion come complete with a fitted cover, hand pump, manual and repair kit.

Cover

Pocket for Ball Pump



Hand Inflation Pump

## StarLock Size Guide

StarLock	5 cm	7.5 cm	10 cm	Seat width*	StarLock	13 cm	Seat width*
26 x 26 cm	ESSL20909-1	ESSL30909-1	ESSL40909-1	22 - 25 cm	26 x 26 cm	ESSL51010-1	22 - 26 cm
26 x 30 cm	ESSL20911-1	ESSL30911-1	ESSL40911-1	22 - 25 cm	26 x 30 cm	ESSL51012-1	22 - 26 cm
26 x 34 cm	ESSL20913-1	ESSL30913-1	ESSL40913-1	22 - 25 cm	26 x 36 cm	ESSL51014-1	22 - 26 cm
30 x 26 cm	ESSL21109-1	ESSL31109-1	ESSL41109-1	26 - 30 cm	30 x 26 cm	ESSL51210-1	27 - 32 cm
30 x 30 cm	ESSL21111-1	ESSL31111-1	ESSL41111-1	26 - 30 cm	30 x 30 cm	ESSL51212-1	27 - 32 cm
30 x 34 cm	ESSL21113-1	ESSL31113-1	ESSL41113-1	26 - 30 cm	30 x 36 cm	ESSL51214-1	27 - 32 cm
30 x 39 cm	ESSL21115-1	ESSL31115-1	ESSL41115-1	26 - 30 cm	30 x 41 cm	ESSL51216-1	27 - 32 cm
34 x 26 cm	ESSL21309-1	ESSL31309-1	ESSL41309-1	31 - 35 cm	36 x 26 cm	ESSL51410-1	33 - 37 cm
34 x 30 cm	ESSL21311-1	ESSL31311-1	ESSL41311-1	31 - 35 cm	36 x 30 cm	ESSL51412-1	33 - 37 cm
34 x 34 cm	ESSL21313-1	ESSL31313-1	ESSL41313-1	31 - 35 cm	36 x 36 cm	ESSL51414-1	33 - 37 cm
34 x 39 cm	ESSL21315-1	ESSL31315-1	ESSL41315-1	31 - 35 cm	36 x 41 cm	ESSL51416-1	33 - 37 cm
34 x 44 cm	ESSL21317-1	ESSL31317-1	ESSL41317-1	31 - 35 cm	36 x 46 cm	ESSL51418-1	33 - 37 cm
39 x 30 cm	ESSL21511-1	ESSL31511-1	ESSL41511-1	36 - 40 cm	41 x 30 cm	ESSL51612-1	38 - 43 cm
39 x 34 cm	ESSL21513-1	ESSL31513-1	ESSL41513-1	36 - 40 cm	41 x 36 cm	ESSL51614-1	38 - 43 cm
39 x 39 cm	ESSL21515-1	ESSL31515-1	ESSL41515-1	36 - 40 cm	41 x 41 cm	ESSL51616-1	38 - 43 cm
39 x 44 cm	ESSL21517-1	ESSL31517-1	ESSL41517-1	36 - 40 cm	41 x 46 cm	ESSL51618-1	38 - 43 cm
39 x 49 cm	ESSL21519-1	ESSL31519-1	ESSL41519-1	36 - 40 cm	41 x 51 cm	ESSL51620-1	38 - 43 cm
39 x 53 cm	ESSL21521-1	ESSL31521-1	ESSL41521-1	36 - 40 cm	41 x 56 cm	ESSL51622-1	38 - 43 cm
44 x 34 cm	ESSL21713-1	ESSL31713-1	ESSL41713-1	40 - 44 cm	46 x 36 cm	ESSL51814-1	43 - 49 cm
44 x 39 cm	ESSL21715-1	ESSL31715-1	ESSL41715-1	40 - 44 cm	46 x 41 cm	ESSL51816-1	43 - 49 cm
44 x 44 cm	ESSL21717-1	ESSL31717-1	ESSL41717-1	40 - 44 cm	46 x 46 cm	ESSL51818-1	43 - 49 cm
44 x 49 cm	ESSL21719-1	ESSL31719-1	ESSL41719-1	40 - 44 cm	46 x 51 cm	ESSL51820-1	43 - 49 cm
44 x 53 cm	ESSL21721-1	ESSL31721-1	ESSL41721-1	40 - 44 cm	46 x 56 cm	ESSL51822-1	43 - 49 cm
49 x 39 cm	ESSL21915-1	ESSL31915-1	ESSL41915-1	44 - 49 cm	51 x 41 cm	ESSL52016-1	48 - 53 cm
49 x 44 cm	ESSL21917-1	ESSL31917-1	ESSL41917-1	44 - 49 cm	51 x 46 cm	ESSL52018-1	48 - 53 cm
49 x 49 cm	ESSL21919-1	ESSL31919-1	ESSL41919-1	44 - 49 cm	51 x 51 cm	ESSL52020-1	48 - 53 cm
49 x 53 cm	ESSL21921-1	ESSL31921-1	ESSL41921-1	44 - 49 cm	51 x 56 cm	ESSL52022-1	48 - 53 cm
53 x 39 cm	ESSL22115-1	ESSL32115-1	ESSL42115-1	49 - 54 cm	56 x 41 cm	ESSL52216-1	53 - 59 cm
53 x 44 cm	ESSL22117-1	ESSL32117-1	ESSL42117-1	49 - 54 cm	56 x 46 cm	ESSL52218-1	53 - 59 cm
53 x 49 cm	ESSL22119-1	ESSL32119-1	ESSL42119-1	49 - 54 cm	56 x 51 cm	ESSL52220-1	53 - 59 cm
53 x 53 cm	ESSL22121-1	ESSL32121-1	ESSL42121-1	49 - 54 cm	56 x 56 cm	ESSL52222-1	53 - 59 cm

\* If placed on a wheelchair, the cushion base edges may fold upwards. Hence the difference between the base measurements and recommended seat width.

## Covers



### Standard Cover


Cover in breathable and flame resistant polyester with 4-way stretch.



### Incontinence Cover

Cover in breathable and flame resistant nylon with 4-way stretch and concealed zipper.

# Lock Pockets – the secret to StarLock

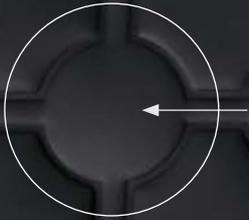

# Setting up StarLock

StarLock can easily be set up for the individual user's need. Once the air level is set, the air can then be locked into the individual cells to create correct support.

## Quick Guide StarLock

**Click or Scan the QR code** to watch the quick guide video for setting the air level for StarLock.

## StarLock Product Page

**Click or Scan the QR code** to find more information about StarLock on StarLock product page, on [etac.com](http://etac.com)





Etac is a world-leading developer of innovative assistive devices and patient handling equipment. Since 1973, we have been committed to improving quality of life for the individual, family members and caregivers.

For the products' intended purpose and latest information, please visit [www.etac.com](http://www.etac.com)

**Etac**



+46 371 587 330



[info@etac.se](mailto:info@etac.se)



[www.etac.com](http://www.etac.com)



Creating Possibilities