

Revolutionized the Industry





The IRIS's adjustable rocker carriage simply aligns the Center of Gravity (CG) over one short base frame

- Maximizes maneuverability, ambulation, and ease of push in tight spaces
- Provides ease of tilt, greater sense of stability for both user and caregiver




Maximum Adjustability, Smart Design

- 8 inches of width adjustment on one frame, without extra parts
- 7 inches of frame depth adjustment on ONE frame with just three seat rails
- Quickly changes with user's needs



Simple, Predictable Factory Set-up, Regardless of Seating

- Takes the guesswork out of aligning the Center of Gravity (CG) with the Center of Rotation (CR)
- Saves time, money, and aggravation for the user, caregiver, clinician, and provider

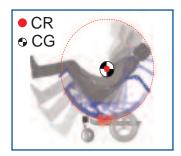

Clinical and Functional Benefits

- Can help reduce startle reflexes by aligning CG with CR
- 60 degrees of tilt, offering better pressure reduction
- Transit tested up to 350 lbs
- Foot activated wheel lock available for caregiver convenience
- 8 back angle adjustments in 5-degree increments
 (-5 to +30 degrees), resulting in seat-to-back angles
 of 85 to 120 degrees in 5-degree increments
- Depth adjustable in 1/2-inch increments, without losing washers
- 4 inches back height adjustment available
- Full range of recline, regardless of tilt position
- Foot propelling set-up available

Degrees of Tilt	-5° to 50°, 0° to 55°, 5° to 60°	
Frame Width Range		14" to 22
Frame Depth Range		15" to 22
Seat to floor height		14.5" to 20.5
Seat to floor height		12.5" to 18.5
Weight Capacity	Std:	250 Lb
	Heavy Duty:	350 Lb
Product Weight*		32 lb:
Back Height	Angle adj. stroller handle	16" to 20
	Standard stroller handle	18", 21", 24
Back Options	Stroller Handle Back Posts (Low	
	Height adjustable with stroller h	nandle, & Reclinin
Tubing Diameter	Seat rail:	1
	Back Post:	1
	Height Adjustable Backrest:	1" and 7/8
Warranty	Frame:	Lifetime
	Accessories:	1 yea
Hanger Options	60°, 70°, 80°, 90°	
	(Swing Away and Heavy Duty),	
	ELR, ALR	
Transit Option	Available on Standard and Hea	vy Duty
Accessories	See Order Form	
Overall Length**		32.5

- * Without: armrest, footrest, wheels and anti-tip
- ${\ensuremath{}^{**}}$ Includes anti-tips and footrests

