


EVA Support Walkers™

HOME

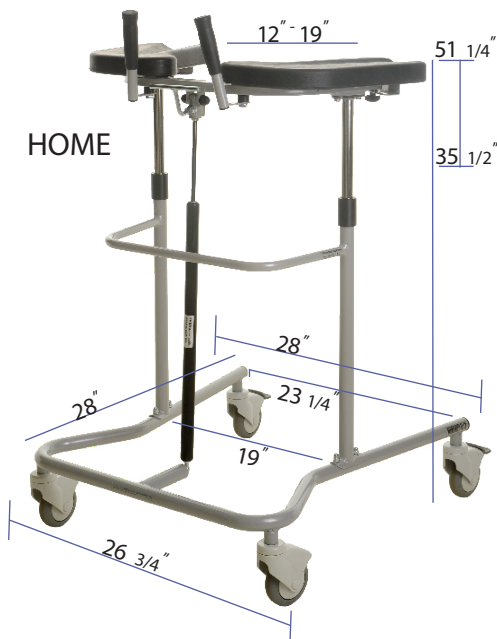

83537
Basket w/Bracket

Sometimes referred to as a “Cardiac Walker”. The Eva® Walker supports the patient’s weight when exercising. Good for Post-op patients when early walking exercise is indicated. Easily adjustable to suit different patient heights. The padded, molded arm pads have a sculpted forearm cutout to keep the arms from sliding off the pads. The arm pads can be removed and are easily cleaned; making it ideal for institutions with a variety of users. Handgrips are a more comfortable round shape and are forward facing, providing a better and more ergonomic grip. Welded steel frame with baked epoxy enamel finish. Non-marring casters with step on wheel locks on the rear casters. Length of adjustable padded armrest is 17 1/2”, and 6 1/2” wide. Width between the armrests can be adjusted from 12”- 19” wide. Step-on locking brakes on two casters.

Adult Pneumatic adjusts from 35 1/2” - 51 1/4” (89/129cm) which should accommodate someone up to 6’5” tall. The easy to release lever mechanism is conveniently located under the front armrest support frame. This unit is ideal when being used with a wide range of patients with it’s easy height adjustment. Includes Standard Handgrips. Weight capacity - 333 lbs.

Directional Casters align the rear casters to track in a straight line to assist the user in controlling the direction of the EVA Walker. Step ON for directional tracking / Step OFF unlocks the directional tracking. NOTE: Caster locks are only on the front casters on models with Directional Casters.

HOME style is for narrow hallways and doors found in the average home.

HOSPITAL style provides a wider entrance for wheelchair use needed by most medical institutions.

- 83502 EVA Pneumatic, Adult, HOSPITAL
- 83507 EVA Pneumatic w/Directional Casters, HOSPITAL
- 83510 EVA Pneumatic w/Hand Brakes, HOSPITAL
- 83540 EVA Pneumatic, SMALL
- 83531 EVA Pneumatic, Adult, HOME
- 83534 EVA Pneumatic w/Directional Casters, HOME
- 83509 EVA Pneumatic w/Hand Brakes, HOME

- 83615 IV Pole w/mounting bracket
- 83617 Portable E & D Oxygen Tank Holder
- 83537 Basket w/Bracket for EVA Walker

EVA ELECTRIC adjusts from 36 1/2” - 50”, is battery operated and will assist the user to a standing position. It has a lifting capacity of (100 kg) 220 lbs. The Pendant Control works independently of the handgrips and is included. Also includes 2 - 12v 2.9 rechargeable batteries for 24V 2.7 Ah hookup. Will operate at 6 minutes per hour continuous use. Min. 50 times up/down between charges. Recommended charging rate is once per week at 12 hours. Batteries and Charger are included. **The 100-240v Charger is CUS, CSA and TO,8A compliant.** Weight capacity 333 lbs. **Push Button Up/Down Handgrips must be purchased separately.**

- 83513 EVA Electric, HOSPITAL
- 83516 EVA Electric, w/ Hand Brakes, HOSPITAL
- 83519 EVA Electric, w/Directional Casters, HOSPITAL
- 83562 EVA Electric, HOME
- 83563 EVA Electric, w/ Directional Casters, HOME
- 83570 EVA Electric, w/Hand Brakes, HOME



83615
IV Pole



HOSPITAL



83617
Oxygen Tank
Holder