

Belly Bag® Urinary Drainage Bag


Discrete Collection



The Belly Bag® is a urine bag that's worn at the waist. It is designed for either males or females who have an indwelling Foley or suprapubic catheter. The bag fastens around the waist by means of a woven belt with a quick release buckle. It has a 1,000 ml capacity and is completely latex-free. An anti reflux valve behind the catheter port prevents reflux urine flow. The Belly Bag® is emptied by means of a drain tube with a twist valve that is in the approximate position of the urethra. Users can easily sleep with the bag in place, thus it can be worn 24 hours a day. In ambulatory patients, the bag is worn under the clothing where a soft backing provides a comfortable fit.



Model is wearing a Belly Bag® filled with 250 ml of water

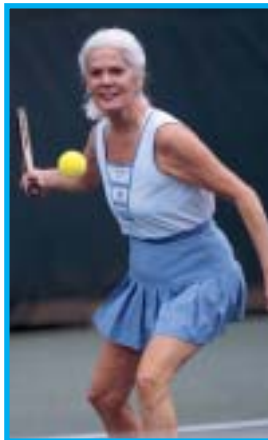

Specifications:

- ◇ 1000 ml capacity
- ◇ Latex-Free
- ◇ Quick release buckle
- ◇ Sterile
- ◇ Anti reflux valve

Benefits:

- ◇ Conceals more effectively
- ◇ Dramatically increases comfort and mobility
- ◇ Eliminates changing between bed and leg bag
- ◇ Helps prevent inadvertent catheter extraction, leg bruises and abrasions

Latex-Free

Belly Bag®		
Product No.	Description	Qty
B1000	1000ml	10
B1000P	1000 ml with sample port	10
B1000CT	1000 ml with 24" drain tube	10
BBDT	Drain tube	10
PO337	Replacement Belt	10

To place an order contact Rüscher Customer Service at **1-800-553-5214**

RÜSCH
A TELEFLEX COMPANY

Teleflex Medical ♦ 2917 Weck Drive
Research Triangle Park, NC ♦ 27709
P: 1-800-553-5214 ♦ F: 770-623-1829
www.rusch.com