


QUICKIE®

SUPERIOR
MANEUVERABILITY

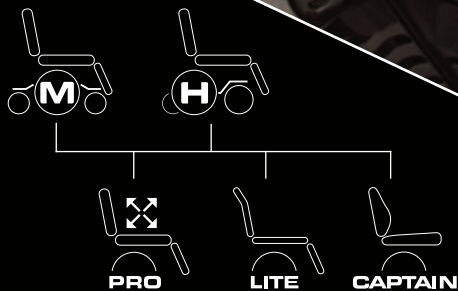

Q500

POWER OF QUICKIE



QUICKIE® delivers the world's most comprehensive range of innovative, technically advanced power chairs, created specifically for your clinical and lifestyle needs.

ADVANCED ELECTRONICS

Controllers are optimized to offer the greatest power for max performance. QUICKIE offers multiple joystick configurations, assignable buttons, on-board programming, and Bluetooth/IR controls. Also available are more battery choices to meet your needs and lifestyle.

SOLID CHASSIS

Chassis is formed from a single sheet metal plate for a stronger, more rigid design. (Others are welded or bolted together.) Steel components undergo automotive undercoat plating to resist rust and corrosion.

OCCUPIED TRANSIT

Tested for ISO-7176 Occupied Transit, included at no charge.

END OF LINE TESTING

After a chair is built, QUICKIE performs a final test for performance, noise, and quality, above and beyond industry requirements.

ENGINEERED FOR CONVENIENCE

QUICKIE engineers are not solely focused on performance and quality, we also thought about convenience and serviceability. QUICKIE Q series power chairs are designed for parts commonality, easy battery access, and simple servicing.

QUICKIE MOTORS

A detailed close-up photograph of a motor's internal gear train. The focus is on a large, polished metal gear in the foreground, with several other gears of varying sizes visible in the background. The lighting is dramatic, highlighting the metallic surfaces and the precision of the engineering.

Our robust motor technology can't be replicated anywhere else. That's because we've spent years working closely with our manufacturer to develop a longer-lasting, highly efficient motor with improved performance, drivability, and range.

DESIGNED TO LAST

Unlike some of our competitors, QUICKIE motors feature metal gears instead of plastic. Polished for low noise and high efficiency, our all-metal gears exceed expectations in the toughest conditions. We go the extra step in polishing those gears because we believe that a little extra effort goes a long way in creating a better ride.

ENGINEERED FOR PERFORMANCE

We've designed our motors to give you the control you need in small spaces without sacrificing the power you want to tackle obstacles. QUICKIE 4-pole motors generate 18% more power than previous generations, plus you'll be able to drive 10% further on the same battery capacity.

MORE EFFECTIVE COOLING

Larger cooling fins enable heat to escape more effectively, reducing your risk of thermal rollback. So you can drive 25% further and faster when combined with the additional boost, even when on an incline or faced with uneven surfaces.

POWER BOOST


Our controller and motor combination work together to provide you with an additional 10 seconds of boost for extra power whenever you need it.

Q500M

MID-WHEEL DRIVE

The QUICKIE Q500 M base balances superb indoor maneuverability with outstanding performance. Its powerful motors combined with all-wheel suspension and QUICKIE's smooth ride, will ensure you have the most comfortable and safe driving experience.

The mid-wheel drive base has the best indoor performance. Its smallest turning radius provides great maneuverability.

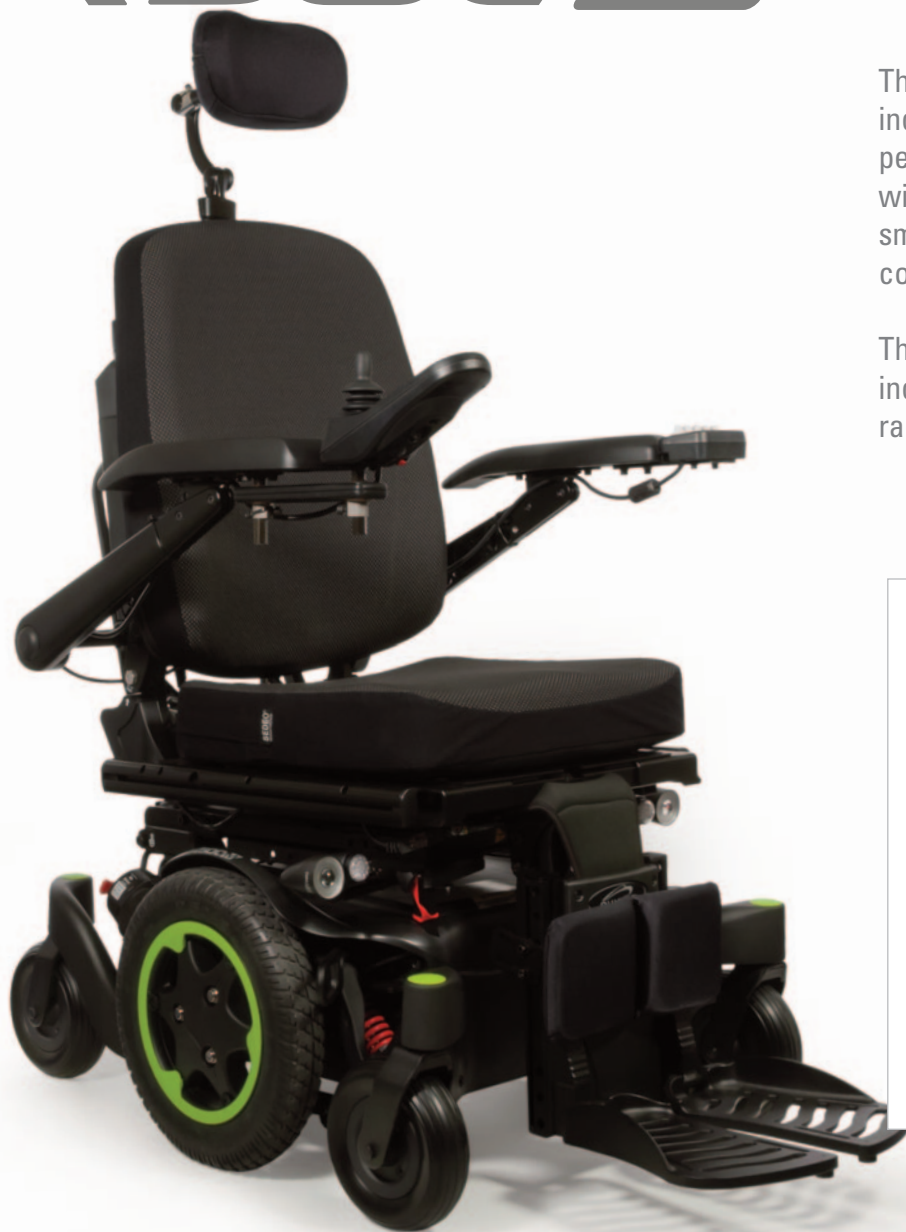



OUTSTANDING PERFORMANCE




ALL-WHEEL SUSPENSION

Enjoy a smooth and stable ride, even over uneven terrain, with QUICKIE's patented all-wheel suspension – six independently suspended wheels give you the smooth ride you're looking for. The Q500 M has a 3" ground clearance and the best traction even when ascending gradients of up to 7°.

GOING THE DISTANCE

With ultra-reliable, high torque 4-pole motors powering speeds of 6 mph and Group 34 batteries standard to maximize your range (up to 14 miles), you can always travel with confidence in the Q500 M.

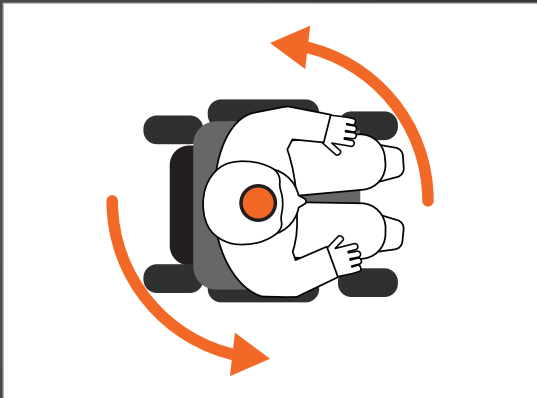

SIMPLE, INTUITIVE DRIVING

Mid-wheel drives offer unmatched maneuverability with a smaller footprint than rear- or front-wheel drives. The Q500 M allows for the most compact turning radius and is extremely intuitive to drive and control.

INDOOR MANEUVERABILITY

The compact dimensions and tight turning radius provide you great access through narrow doorways of as little as 25", on lifts and around confined spaces. Easily transfer in and out and effortlessly access vehicles and under tables with the lower than ever before seat-to-floor height (from 16.5" with powered center of gravity tilt).

Q500H

HYBRID-WHEEL DRIVE

QUICKIE Q500 H brings you a unique blend of rear-wheel drive and mid-wheel drive technology delivering impressive outdoor performance coupled with indoor maneuverability. Bringing you the best of both worlds and exceptional comfort.





A NEW GENERATION OF POWER



STABLE RIDE

Intelligent, active rear anti-tips ensure all 6 wheels remain in contact with the ground for enhanced traction and stability. The Q500 H's fixed front casters prevent the chair from tipping forward while the dynamic suspension system keeps both drive wheels and casters in contact with the ground.

INCREASED MANEUVERABILITY

Unlike conventional rear wheel drive bases the Q500 H's drive wheel position ensures a smaller turning radius for increased indoor maneuverability but without compromising stability and outdoor performance.



EXCEPTIONAL RIDE COMFORT

The innovative, dynamic suspension system actively absorbs impacts over thresholds and rough terrain to enhance ride comfort and seat stability and to reduce fatigue and ultimately extend the time a user can sit in the chair.

PERFORMANCE

The large front casters not only provide a more stable ride, but they allow greater ability to climb obstacles. Additionally, these large casters absorb more impact giving you both performance and a smoother ride.

CHOOSE YOUR SEAT FRAME

SEDEO® PRO™

Truly modular in design, the SEDEO PRO SEAT FRAME provides a new level of functionality, comfort and design. For moderate to complex clinical needs, the PRO combines a highly adjustable seat frame with an excellent choice of postural support options.



SEDEO® LITE™

The SEDEO LITE is an adjustable seat frame with a sleek look. Its minimal design alleviates weight capacity concerns when using aftermarket equipment. Traditional back canes are width-adjustable from 15" to 20" and easily fit JAY® backs. An array of options cover mild to moderate clinical needs, and intuitive adjustability meets individual requirements.



ELEVATE AND DRIVE WITH C-ME®

No more waiting around or driving at a reduced speed when elevated. C-Me enables you to elevate up to 12" while traveling up to 3 mph! You can lift or descend at the touch of a button with no need to stop while transitioning. Just rise-up and go at the same time!



CAPTAIN’S SEAT

The Captain’s Seat provides a relaxed sitting position and features a semi-recline, fold-down back, adjustable seat depth, and convertible seat pan that allows you to use your preferred cushion. Available with a 50° tilt and power centermount legrest.



SEAT FRAME COMPARISON	SEDEO PRO	SEDEO LITE	CAPTAIN’S SEAT
Clinical Need	Med to High	Med to High	Low
CryoFluid™ (Available only on JAY® Fusion & Syncra Seating)	●	●	●
Powered Recline with 4" of Shear Reduction (150°)	●	X	X
Semi Power Recline Option (130°)	●	X	X
Manual Recline Option	●	X	●
Power 50° Center of Gravity Tilt	●	●	●
Power Legrest	●	●	●
12" Elevate	●	●	●
Seat Width	16" to 22"	15" to 20"	16", 18", 20"
Seat Depth	14" to 22"	16" to 22"	16" to 22"
Back Height	21" to 30"	18" to 30"	25"

SEDEO® PRO™

HIGHLY CUSTOMIZABLE SEATING

The SEDEO PRO is a fully adjustable and configurable seating frame with a wide range of comfort cushions and positioning options. With its diverse selection of armrests, armpads, legrests, seat and back cushions, headrests, and positioning accessories, SEDEO PRO offers great flexibility to adapt to individual needs.

QUICK AND EASY TO ADJUST

We've made it easy for healthcare professionals to create the right seating posture and maintain it for changing needs.

- Intuitive adjustments using just one tool
- A wide range of adjustments for varying and changing needs
- Seat width and depth, back height, and armrest width all can be adjusted and repositioned to create the best support in the shortest amount of time



WIDE ARRAY OF OPTIONS



TILT, ELEVATE, AND RECLINE

The SEDEO PRO offers 50° of tilt with automatic COG adjustment for improved weight distribution and stability. In addition to the low seat-to-floor height, enjoy 12 inches of elevate, rated at up to 300 lbs!



BACK SUPPORT OPTIONS

Choose from 4 options: (1) fixed, (2) user-adjustable manual semi-recline, (3) 130° powered semi-recline or (4) a 150° powered recline back. The 150° powered recline back offers 4" of anti-shear to minimize the sliding between you and the cushion.

LEG AND FOOT POSITIONING

The SEDEO PRO offers a full array of foot support options. You can choose a power elevating and articulating legrest which allows for optimal positioning and comfort. Power elevating legrests are available as centermount or individual swing-away versions. Alternatively, you can choose a manual elevating swing-away legrest or traditional fixed angle style in both swing-away or centermount.

EXCEPTIONAL, VERSATILE POSITIONING

Many positioning elements are available such as headrests, lower extremity supports, swing-away lateral thoracic supports, thigh supports, or knee supports. All positioning elements are adjustable and detachable or can be swung out of the way for easier transfers.

FLIP-BACK ARMRESTS

The easy to operate flip-back armrest gives you unencumbered transfer space. Choose from a removable or reclining version. A single post removable armrest is also available. All are compatible with four armpad options.



SEDEO® LITE™

SIMPLE, ADJUSTABLE SEATING

Sometimes less is more. The SEDEO LITE is an adjustable seat frame giving you a simple, sleek look. Its minimal design alleviates weight capacity concerns when using aftermarket equipment. As one of the lightest seating systems in the market, the SEDEO LITE has a starting weight of 266 lbs, making it a great solution for many hitch mounted lifts.



EASY TO ADJUST

Effortlessly make a wide range of adjustments using just one tool!

- To accommodate changing requirements, seat widths can adjust between 15" and 20", and seat depth can adjust between 16" and 22".
- The easy-to-operate, flip-back armrest gives you unencumbered transfer space.
- With 5" of armrest height adjustability and a large array of adjustable options available on the SEDEO LITE, it is easy find the right fit.
- The side seat rail allows for additional mounting solutions such as transfer bars, thigh supports and knee supports with infinite adjustment.
- The back canes are width-adjustable from 15" to 20".



FIND THE PERFECT FIT WITH NOTHING EXTRA TO WEIGH YOU DOWN



TILT AND ELEVATE

The SEDEO LITE offers 50° of tilt with automatic COG adjustment for improved stability. In addition to the low seat-to-floor height, enjoy 12 inches of elevate, rated up to 300 lbs!



SEATING AND SUPPORT OPTIONS

- The traditional back canes easily fit a range of JAY® backs, and the clamp hardware allows infinite height adjustment.
- Alternatively, you can choose a height-adjustable sling back with adjustable straps.
- Choose from a selection of Whitmyer® or SEDEO® headrests – compatible whether you opt for a traditional back mount or the integrated mount with the sling.



FOLD-DOWN BACK

With transportation in mind, the lightweight seat is complimented with a fold-down or removable back to provide one of the most compact power seating systems in the market.

LEG AND FOOT POSITIONING

The SEDEO LITE offers a full array of legrest options. Choose from multiple centermount options such as fixed, power, or angle-adjustable.



JAY® SYNCRA™ SEATING



POSITIONED FOR RECLINE

JAY® Syncra™ is a cohesive and attractive seating system specifically designed for performance with tilt and recline systems. Syncra has a width and height adjustable back, a cushion that fits mild to complex seating needs, and positioning accessories for support and stability.

JAY® Syncra™ is the optimal seating for power wheelchair users who have a history of pressure injuries.

DESIGNED EXCLUSIVELY FOR
THE QUICKIE® SEDEO® PRO™

JAY CRYOFLUID™

The JAY Syncra cushion is also available with the revolutionary, patented JAY Cryo™ Technology. JAY Fluid evenly distributes pressure across the seating surface and around bony prominences while the **CryoFluid™ actively cools** the skin for up to eight hours.* The combined technology offers a dramatic improvement in tissue injury prevention.



*Internal testing data at 25°C. Results may vary.

JAY® & WHITMYER® SEATING OPTIONS

A wide range of JAY cushions and backs are optimized to fit the SEDEO PRO or SEDEO LITE for advanced positioning capabilities. All the style, comfort and adjustability you've come to love with WHITMYER has been carried over to SEDEO Seating, with a headrest post design that integrates seamlessly into the back support.

Designed for pelvic stability, skin protection and comfort, JAY cushions and backs are low maintenance, highly personalized and available in a wide variety to meet your needs.



JAY Union®



JAY J3® Back



JAY Fusion®

WHITMYER® HEADRESTS

Purposefully engineered to address a range of needs from basic support to multiple-point control, WHITMYER headrests have five mounting options and a variety of pad shapes and sizes for superior head positioning and comfort.



Adjust-A-PLUSH w/LINX



Contoured Cradle w/LINX²



PLUSH w/AXYS



PLUSH w/Cobra

R-NET ADVANCED JOYSTICK

The Q500 is available with the superior R-net Advanced Joystick. A large, easy-to-read display – largest in the market – and simple to control paddle switches allow for intuitive operation.

CONNECT


Seamlessly take control of your environment using the integrated infrared functionality. Use your joystick like a remote control for your TV and select home entertainment systems.

BLUETOOTH® WIRELESS TECHNOLOGY

Control your phone, computer, tablet and all Bluetooth® devices with this useful option. Connect and control up to four devices. Compatible with Android™ and Apple® iOS devices.

QUICKIE PATENTED ASSIGNABLE BUTTONS

Change the function of any button on your joystick to any command you desire, thanks to QUICKIE assignable buttons. You can also assign a second (long press) function to buttons, giving you quicker, simpler and smoother access to functions you commonly use without having to use your joystick.

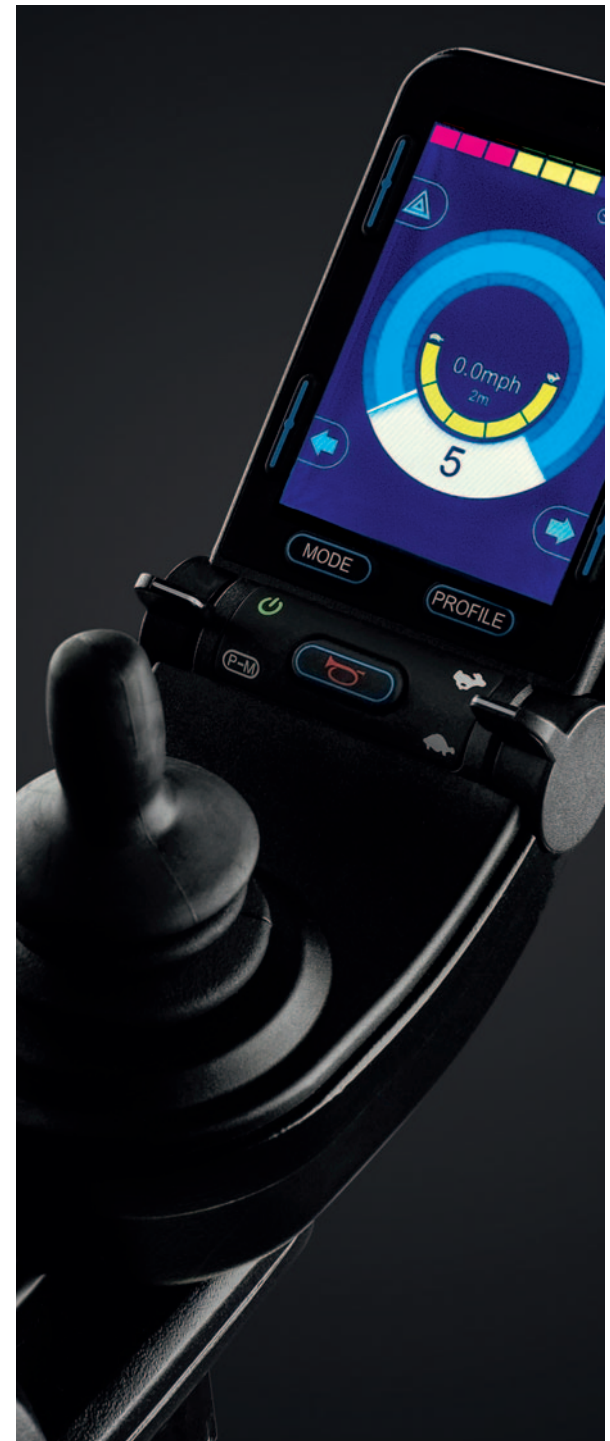

SWITCH-IT™ INNOVATIVE CONTROLS

CTRL+5®

Conveniently access five seating functions with the Ctrl+5. It allows you to go straight to the function you need at the touch of a button.

Available in push-button (shown) or toggle.



ALTERNATIVE DRIVE CONTROLS

If driving with a traditional joystick is not the best option, alternative drive controls allow you to fully experience your chair's driving capabilities and benefits.

MICROGUIDE— Light-touch proportional control.

If you could benefit from light force requirements but need the feedback of joystick movement, the MicroGuide is the ideal choice.

VERSAGUIDE— Compact proportional control.

Fully proportional, durable, and compact, the VersaGuide can be mounted virtually anywhere. It comes with two on/off, mode SmartJacks which automatically adapt to both proximity and mechanical switches.

DUAL PRO®— Truly proportional head array for intuitive driving.

The Dual Pro is a truly proportional head array with fine-tuned adjustments and customization. The Dual Pro can be programmed conveniently from the back of the head array and customized to your strengths and capabilities.



LINK-IT™ MOUNTING SYSTEM

Need to use additional switches or specialty controls? No problem! Easily connect them using the universal Link-It Mounting System.

Using a combination of available modules and rods, create endless combinations. It's easy to position and ensures your wheelchair maintains its sleek look. Whether it's an extra controller on the armrest, a push button on the back support or a switch on the headrest, just link-it.



CAPTAIN SEAT RELAXED COMFORT



The Captain's Seat provides a relaxed sitting position and features a semi-recline, fold-down back, adjustable seat size, and convertible seat pan that allows you to use your preferred cushion. Available with a 50° tilt and power centermount legrest.

- Depth adjustable seat size, legrest placement, and armrest width
- Convertible seat pan for flexibility in cushion selection
- Semi-recline and fold-down back
- Optional oxygen, crutch, and walker holders




DISTINGUISHED STYLE

With the Q500 series, we've taken style to the next level. With a matte black chassis finish, hidden wiring and premium upholstery, the Q500 has a sleek and contemporary design.

The Q500 comes with a Matte Black shroud and your choice of six contemporary color inlays to beautifully accent the matte black bodywork. Fancy a change? Simply replace the color inlays. Styling your powered wheelchair has never been so simple.




Additionally, the Q500 M offers your choice of eight different shroud colors to reflect your personal style. (Color accents are only available with the Matte Black shroud.)

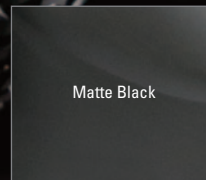


Matte Black



Candy Apple Red

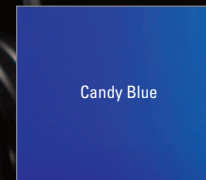


Candy Blue

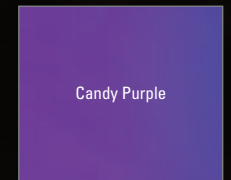

Candy Purple



Apple Green

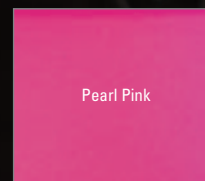

Pearl Pink



Silver



American Flag

IMPORTANT NOTE— Print processes cannot accurately represent painted colors. Actual colors may vary. For more accurate representations of colors, refer to the QUICKIE/ZIPPIE Color Guide or a QUICKIE Color Ring.

	Q500 M	Q500 H		
HCPCS codes	Group 3 – Standard Duty No Power Seating: K0848 Single Power Seating: K0856 Multi Power Seating: K0861	Group 3 – Standard Duty No Power Seating: K0848 Single Power Seating: K0856 Multi Power Seating: K0861	Group 3 – Heavy Duty No Power Seating: K0850 Single Power Seating: K0858 Multi Power Seating: K0862	Group 4 – Standard Duty No Power Seating: K0868 Single Power Seating: K0877 Multi Power Seating: K0884

	Q500 M	Q500 H
Max. user weight	300 lbs.	300 lbs. / 350 lbs. HD
Seat width	15" to 22"	16" to 22"
Seat depth	14" to 22"	14" to 22"
Seat-to-floor height	16.5" to 19"	16.5" to 19"
Back height	18" to 30"	21" to 30"
Total width	24"	24.5"
Total length (without footplates)	34.5"	36"
Min. turning radius (without footplate)	18.5"	25.5"
Max. speed	6 mph	6 mph (standard) / 5 mph (HD) / 8 mph (optional)
Battery size	Group 22 / Group 34	Group 24 / Group 34
Overall weight	No Power Seating: starting at 266 lbs. Single Power Seating: starting at 291 lbs. Multi Power Seating: starting at 338 lbs.	No Power Seating: starting at 343 lbs. Single Power Seating: starting at 367 lbs. Multi Power Seating: starting at 381 lbs.
Drive wheel position	Mid	Mid/Rear (Hybrid)
Max. Ground clearance	3"	3"
Wheel options	13" drive wheels, 7" casters	13" or 14" drive wheels, 9" casters
Max. range	Up to 14 miles	Up to 17.7 miles (w/ GP24 batteries and 5 mph motors)
Transit Securement Brackets	ISO 7176-19:2008 Transit Approved	ISO 7176-19:2008 Transit Approved

See order form as specifications may vary based upon configuration.



Sunrise Medical (US) LLC, Fresno, CA 93727
MK-129555 Rev. E ©9.2020

800.333.4000 www.SunriseMedical.com

