



QUICKIE®

THE ULTIMATE  
EXPERIENCE

QUICKIE  
Q700M

# THE ULTIMATE EXPERIENCE

QUICKIE® delivers the world's most comprehensive range of innovative, technically advanced power chairs, created specifically for your clinical and lifestyle needs.



Resulting from years of experience, the QUICKIE® Q700™ M offers the highest performance and technology, providing you with The Ultimate Experience!

## HIGH-PERFORMANCE BASE

Advanced mid-wheel drive technology gets you where you want to go. SpiderTrac® 2.0 suspension and our patented anti-pitch system provide excellent outdoor performance and smooth curb climbing.

## TWO INTELLIGENT SEATING OPTIONS

### SEDEO® ERGO™

A highly specialized and technical seating system that allows for precision fitting for complex clinical needs. The SEDEO ERGO includes top-of-the-line features, such as Biometric Repositioning and programmable seating.

### SEDEO® PRO™

An extremely configurable and versatile seating system, the SEDEO PRO can be easily customized to suit individual requirements and quickly adapted to changing needs. It features a wide variety of options to fit almost everyone's requirements.

## INNOVATIVE CONTROLS

Seamlessly connect with your environment, your smart phone, computer and other electronics with a choice of industry-leading, intuitive and easy-to-operate controls, including infrared and BLUETOOTH® technology.

## DISTINGUISHED STYLE

The Q700 M has a streamlined, minimalistic design contrasted with sleek, black bodywork and striking color inlays. Every detail from the upholstery down to the wheels has been designed to turn heads.







## MID-WHEEL DRIVE PERFORMANCE, EVOLVED

The Q700 M Mid-Wheel Drive base sets the bar for powered wheelchair technology, inside and out. Travel further with its superior drive performance. Feel secure over challenging terrain with SpiderTrac 2.0 suspension and its independent movement of all 6 wheels. It's a driving experience like never before.

## TRAVEL WITH CONFIDENCE

Ultra-reliable, high torque 4-pole motors power speeds of 5, 6 and 8 mph. A choice of Group 24 or Group 34 batteries maximize your range. With occupied transit securement brackets included at no charge, you can always travel with confidence in the Q700 M.




# OUTDOOR PERFORMANCE



Climb curbs of up to 3" and enjoy a smoother, more responsive ride with the revolutionary SpiderTrac 2.0 suspension. With its independent movement of all six wheels, the SpiderTrac suspension gives the Q700 M the ability to safely crawl up and down steep obstacles. The performance springs allow for more wheels in contact with the ground at any time – even when ascending gradients of up to 9°. *Limitations may apply.*

*Please read the owner's manual for safe operating instructions.*

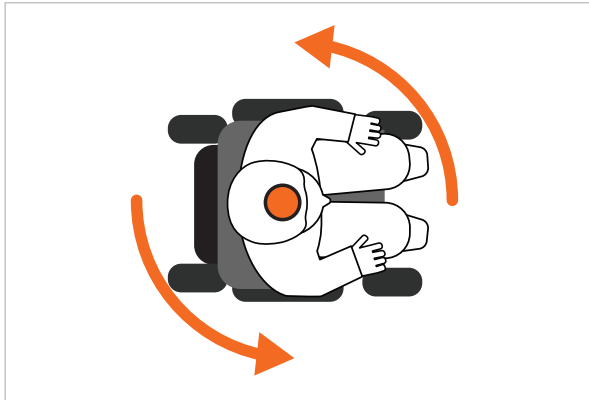


# INDOOR MANEUVERABILITY

Turning on its own axis, enjoy intuitive driving and excellent maneuverability in tight compact spaces with this narrow 25" base and a 23.5" turning radius. Also, easily roll under tables and access Wheelchair Adaptive Vehicles (WAVs) with a low minimum seat height of just 16.5".



WITH SURETRAC®

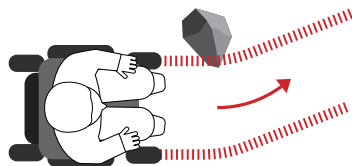

WITHOUT SURETRAC®

# STAY ON TRACK WITH SURETRAC®

The optional SureTrac system automatically detects and corrects the smallest deviation from your input direction. By providing a more accurate drive, it will reduce the number of physical or cognitive actions required to keep you on track, making for a less demanding drive (especially if using specialty controls).

# CHOOSE YOUR SEATING

## SEDEO® ERGO™

Designed for comfort and aesthetics, the SEDEO ERGO covers the most demanding clinical needs. With Biometric Repositioning and a full range of programmable seating positions, it is the best ergonomic seating from QUICKIE.



## SEDEO® PRO™

Truly modular in design, the SEDEO PRO seating system provides a new level of functionality, comfort and design. For moderate to complex clinical needs, the PRO combines a highly adjustable seat frame with an excellent choice of postural support options.



## ELEVATE AND DRIVE WITH C-ME® AVAILABLE ON BOTH ERGO AND PRO

No more waiting around or driving at reduced speed when elevated. C-Me enables you to elevate up to 12" while traveling up to 3 mph! (Weight limitations may apply.\*) You can lift or descend at the touch of a button with no need to stop while transitioning. Just rise-up and go at the same time!





## SEATING SYSTEM COMPARISON

	<b>SEDEO ERGO</b>	<b>SEDEO PRO</b>
<b>Clinical need</b>	High	Med to High
<b>Seat width</b>	16" to 22"	16" to 22"
<b>Seat depth</b>	16" to 22"	16" to 22"
<b>Backrest height</b>	21" to 28"	21" to 30"
<b>Biometric Repositioning</b>	●	
<b>Memory seating</b>	●	
<b>Powered recline</b>	165°	150° / 130°
<b>Anti-shear backrest</b>	7"	4" (optional)
<b>Pressure relief alert reminders (PRAR)</b>	●	
<b>Tilt angle</b>	50°	50°
<b>Lift</b>	12"	12"
<b>JAY &amp; WHITMYER</b>	●	●
<b>Remote Seating App</b>	●	



**QUICKIE Q700 M**  
**SEDEO ERGO SEATING SYSTEM**

# SEDEO® ERGO™ INTELLIGENT SEATING

It's the first-class seat that's designed specifically to fit you. The SEDEO ERGO with patented Biometric Repositioning is more effective than any other seating system when it comes to independent repositioning throughout the day.

## UPHOLSTERY OPTIONS

- Black/Grey Leatherette for a premium look
- 3DX™ Spacer Fabric for maximum breathability
- 3DX™ Spacer Fabric with Lycra® outer layer for lowest shear

## EASY TRANSFER

The seat cushion is divided into 3 elements. The two thigh bolsters can be swung out and lifted off, ensuring a comfortable transfer. The center seat cushion can also be substituted with a JAY® cushion.

## MEMORY POSITIONS

With memory seating, you can save up to six programmable positions and transition with the simple press of a button.







## ADJUSTABLE TO FIT YOU

The premium headrest (1) can be easily adjusted or removed. All seat widths and depths (2) are infinitely adjustable with just one tool. Also, the seat and backrest width (3) can be adjusted independently of each other. Even the armrest height and angle (4) are simple to adjust.

### PRESSURE RELIEF REMINDERS

Pressure relief reminders tell you when it's time to change position. A log can also be kept for clinicians to review.

### PROGRAMMABLE MEMORY SEATING

The SEDEO ERGO has 6 uniquely customizable seating positions that can be programmed to any position you desire. Whether you want to lie down, cook, or watch TV – every position can be completely customized to you.

### FOOTPLATE OPTIONS

Choose from 2 sizes of footplates to meet your unique needs, and select the optional ability to lower the footplate all the way to the ground, making transfers easier.



## BIOMETRIC REPOSITIONING

Moving out of a reclined position often requires use of the tilt function to regain your proper seated posture. With the patented Biometric Repositioning, when you move from reclined to seated, the SEDEO ERGO will first pass through an intermediary tilted position, using gravity to realign your pelvis in one smooth, effortless motion.



## ANTI-SHEAR BACKREST

Transitioning in and out of recline doesn't need to be uncomfortable. The SEDEO ERGO's power recline is supported by 7" of anti-shear to better maintain the positioning of the headrest, and even chin control, during recline.

# SEDEO® PRO™ HIGHLY CUSTOMIZABLE SEATING

Fully adjustable and configurable seating frame with a wide range of cushions and positioning options. With its diverse selection of armrests, armpads, leg rests, seat and back cushions, headrests and positioning accessories, SEDEO PRO offers great flexibility to adapt to individual user needs.



## FLIP-BACK ARMRESTS

The easy to operate flip-back armrest gives you unencumbered transfer space. Choose from a single-post or reclining version. All are compatible with four armpad options.



## QUICK AND EASY TO ADJUST


Creating the right seating posture is easy to do for the healthcare professional.

- Intuitive adjustments using just one tool.
- A wide range of adjustments can be made for your varying and changing needs.
- Seat width and depth, back height and armrest width all can be adjusted and repositioned to create the best support in the shortest amount of time.

## WIDE ARRAY OF OPTIONS

### CHOICE OF BACKREST OPTIONS

A variety of backrest options are available for the SEDEO PRO. Choose from 4 options: (1) fixed, (2) user-adjustable manual semi-recline, (3) powered semi-recline or (4) a full powered recline backrest. The full powered recline backrest offers 4" of anti-shear, making sure to minimize the sliding between you and the cushion.

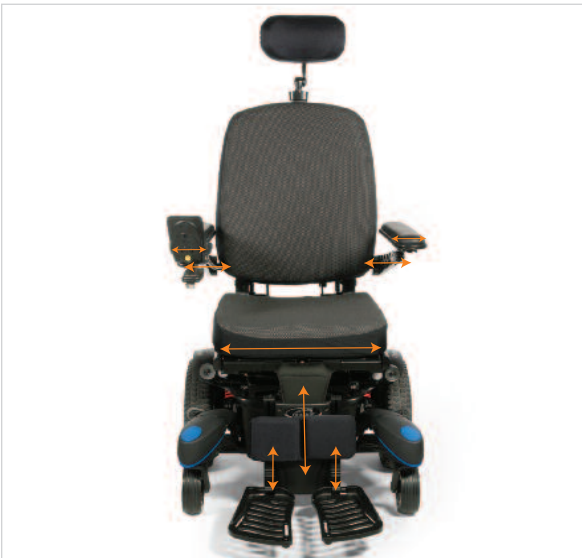

### LEG AND FOOT POSITIONING OPTIONS

The SEDEO PRO offers a full array of footrest options. You can choose a power elevating and articulating leg rest which allows for optimal positioning and comfort. Power elevating leg rests are available as center mount or individual swing-away versions. Alternatively, you can also choose a manual elevating swing-away leg rest or traditional fixed angle style in both swing-away or center mount.

### EXCEPTIONAL, VERSATILE POSITIONING


Many positioning elements are available such as headrests, lower extremity supports, swing-away lateral thoracic supports, thigh supports or knee supports. All positioning elements are adjustable for prime positioning and are detachable or can be swung out of the way for easier transfers.

### JAY® AND WHITMYER® COMPATIBLE

A wide range of JAY cushions and backs were optimized to fit the SEDEO PRO for advanced positioning capabilities. All the style, comfort and adjustability you've come to know and love with WHITMYER has been carried over to the SEDEO PRO with a new design that integrates the headrest post seamlessly into the PRO back. Several options are shown on page 17.





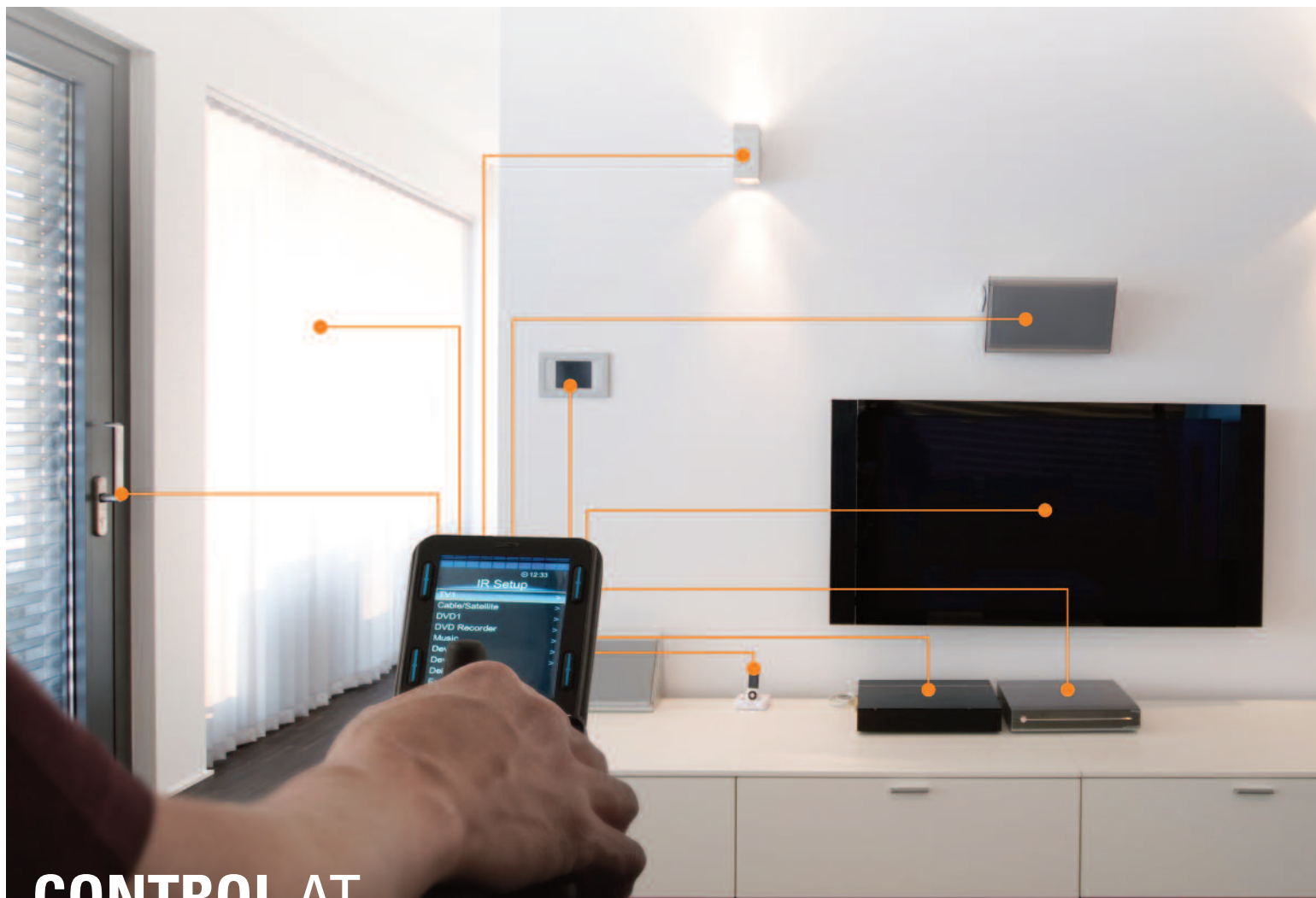

# CONTROL AT YOUR FINGERTIPS

Be the master of your environment with the Q700 M Series and its advanced controller technology. Not only can you control your Q700 M the way you want, but also your phone, computer, tablet and all Bluetooth® devices, including Android™ and Apple® iOS devices.

# R-NET ADVANCED JOYSTICK

Every Q700 M Series comes with the superior R-net Advanced Joystick. A large, easy-to-read display – largest in the market – and simple to control paddle switches allow for intuitive operation.

## CONNECT


Seamlessly take control of your environment using the integrated infrared functionality. Use your joystick like a remote control for your TV and select home entertainment systems.

## BLUETOOTH® WIRELESS TECHNOLOGY

Control your phone, computer, tablet and all Bluetooth® devices with this useful option. Connect and control up to four devices. Compatible with Android™ and Apple® iOS devices.

## QUICKIE PATENTED ASSIGNABLE BUTTONS

Change the function of any button on your joystick to any command you desire, thanks to QUICKIE assignable buttons. You can also assign a second (long press) function to buttons, giving you quicker, simpler and smoother access to functions you commonly use without having to use your joystick.

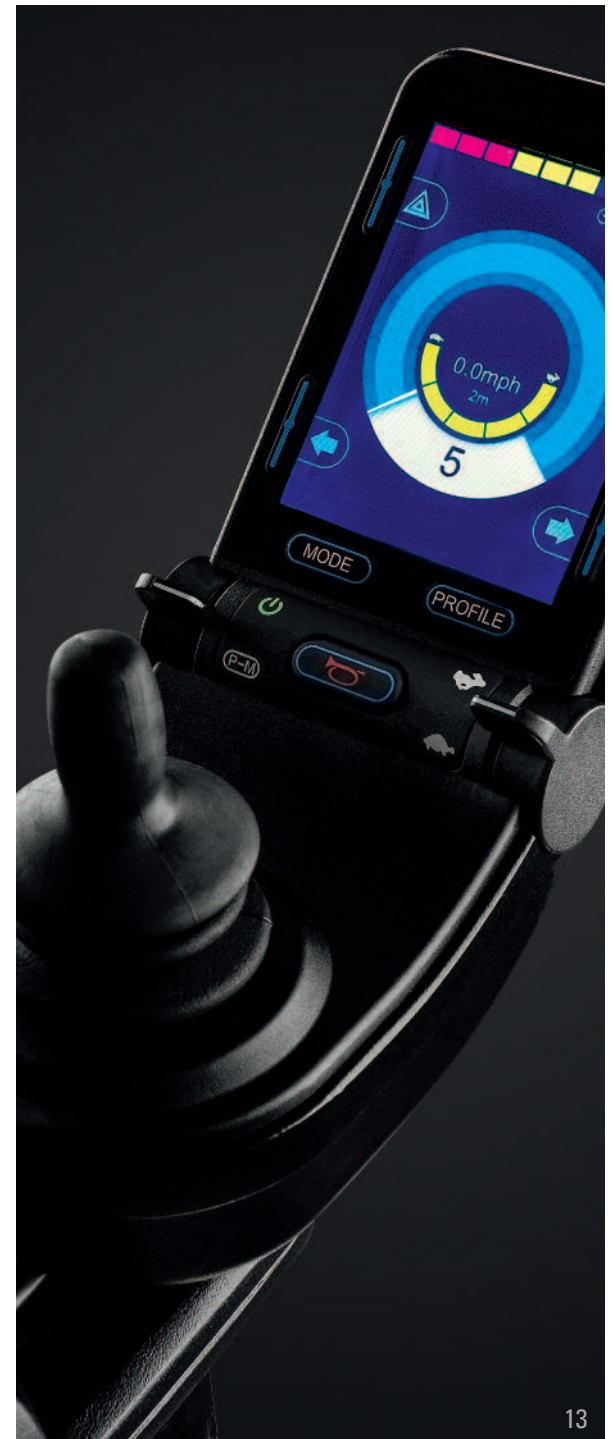

# SWITCH-IT™

## INNOVATIVE CONTROLS

### CTRL+5®

Conveniently access five seating functions with the Ctrl+5 – a perfect combination with Memory Seating. It allows you to go straight to the function you need at the touch of a button.

Available in push-button (shown) or toggle.



### ALTERNATIVE DRIVE CONTROLS

If driving with a traditional joystick is not the best option, alternative drive controls allow you to fully experience your chair's driving capabilities and benefits.

#### MICROGUIDE

**Light-touch proportional control.** Proportional input devices offer the fine control you need for independent access of your world. If you could benefit from light force requirements but need the feedback of joystick movement, the MicroGuide is the ideal choice.



#### VERSAGUIDE

**Compact proportional control.** The VersaGuide is a fully proportional and durable compact joystick that can be mounted virtually anywhere. Its unlimited mounting options are what give the versatile VersaGuide its name. The VersaGuide comes with two built-in jacks for mode and on/off, called SmartJacks, which are unique to the VersaGuide. The SmartJack sensor input will automatically adapt to both proximity and mechanical switches.



#### DUAL PRO®

**Truly proportional head array for intuitive driving.** The Dual Pro is a truly proportional head array with the ability to fine-tune adjustments and customization to provide you an easy-to-learn, precise, and intuitive driving experience. The Dual Pro offers a full range of acceleration and deceleration, 360° turning range and veering control for superior drivability. The Dual Pro can be programmed conveniently from the back pad of the head array and customized to your strengths and capabilities.







# SWITCH-IT™ REMOTE SEATING APP

AVAILABLE EXCLUSIVELY ON SEDEO ERGO SEATING

## TRACK. ACT. THRIVE.


The Remote Seating app plays an important role in pressure management. With the app, you can self-monitor your pressure relief progress and securely share the data with your caregiver network. The app tracks key seating positions such as tilt, recline and leg elevation, as well as the duration of these positions.

## DATA. WHEN YOU WANT IT. WHERE YOU WANT IT.

With the one-touch app, there's no need to log into a cumbersome portal, just connect to your chair. The data stays with your chair and is easily accessible via your phone or tablet. You can also email reports instantly to your choice of clinicians and doctors.

- Secure data
- Easy access
- Bluetooth® connection
- Live progress reports
- View data directly on your device
- Assists in pressure relief plans
- Pressure relief alert reminders
- Compatible with Android™ and Apple® iOS

*(Refer to owner's manual for compatible versions.)*

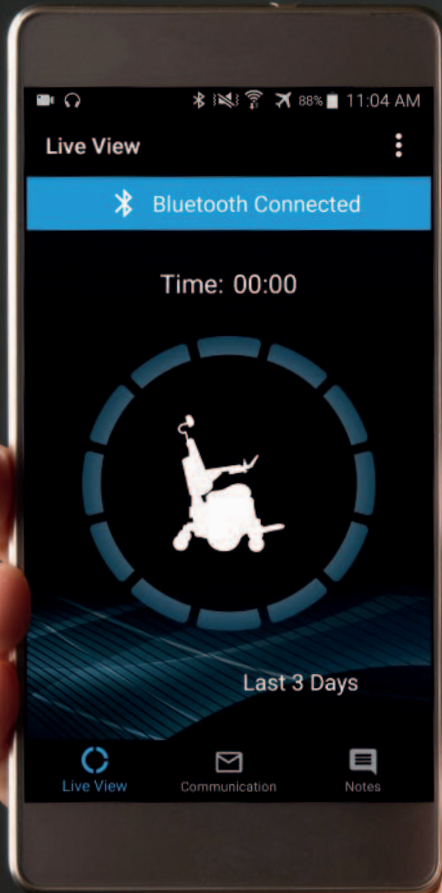

## TAKE CHARGE OF YOUR SKIN HEALTH.

# LINK-IT™ MOUNTING SYSTEM

Need to use additional switches or specialty controls? No problem!  
Easily connect them using the universal Link-It Mounting System.

Using a combination of available modules and rods, create endless combinations.  
It's easy to position and ensures your wheelchair maintains its sleek look.  
Whether it's an extra controller on the armrest, a push button on the backrest  
or a switch on the headrest, just Link-It.



Knee Mount



# JAY® & WHITMYER® SEATING OPTIONS

A wide range of JAY cushions and backs and WHITMYER headrests are available with the SEDEO ERGO and PRO seating systems.

## JAY® CUSHIONS & BACKS

Designed for pelvic stability, skin protection and comfort, JAY cushions and backs are low maintenance, highly customizable and available in a wide variety to meet your needs.



JAY Union®



JAY J3® Back



JAY Fusion®

## WHITMYER® HEADRESTS

Purposefully engineered to address a range of needs from basic support to multiple-point control, WHITMYER headrests have five mounting options and a variety of pad shapes and sizes for superior head positioning and comfort.



Adjust-A-PLUSH w/LINX



Contoured Cradle w/LINX²



PLUSH w/AXYS



PLUSH w/Cobra







# DISTINGUISHED STYLE

With the Q700 M series, we've taken style to the next level. With a matte black chassis finish, hidden wiring and premium upholstery, the Q700 M has a sleek and contemporary design.

The Q700 M also comes with a choice of six contemporary color inlays to perfectly complement its striking matte black bodywork. Fancy a change? Simply replace the color inlays. Styling your powered wheelchair has never been so simple.

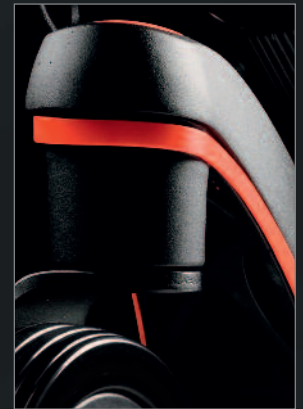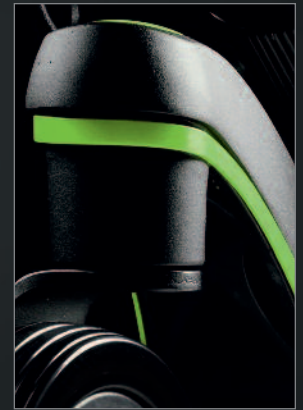

# Q700 M

See order form as specifications may vary based upon configuration

HCPC codes	No Power Seating Options	Single Power Seating Options	Multi Power Seating Options
Group 3 Standard Duty	K0848	K0856	K0861
Group 3 Heavy Duty	K0850	K0858	K0862
Group 4 Standard Duty	K0868	K0877	K0884
Max. user weight	300 lbs* (350 lbs for HD)		
Seat width	16" to 22"		
Seat depth	16" to 22" <i>Note: Back cushion can cause 2" loss of effective depth.</i>		
Seat-to-floor height	16.5" to 19"		
Backrest height	21" to 30"		
Total width	25"		
Total length	36.1"		
Turning radius (includes footplate)	23.5"		
Seat memory functions (ERGO only)	Seated, Pressure Relief, Supine, Feet Above Heart, Transfer, C-Me		
Max. speed	6 mph (standard), 5 mph (optional), 8 mph (optional)		
Battery size	Group 34 or Group 24		
Overall weight	Starting at: 325 lbs (w/o batteries), 410 lbs (w/batteries)		
Drive wheel position	Mid		
Ground clearance	3"		
Wheel options	14" drive wheels; 6" casters		
Max. range	Up to 18 miles (w/Group 24 batteries and 6 mph motors)		
Transit Securement Brackets	ISO 7176-19:2008 Transit Approved		
Warranty	Lifetime Limited		

\*C-Me at 3 mph has a max. user weight of 250 lbs.



Sunrise Medical (US) LLC, Fresno, CA 93727  
MK-129475 Rev.C ©10.2018

800.333.4000 [www.SunriseMedical.com](http://www.SunriseMedical.com)

