

AQUATEC®

Ergo Family



Modular personal care solutions

SAFE | INDEPENDENT | ERGONOMIC



Customer insight-led design to enhance independence, safety and flexibility.

Personal care is highly private activity and is an important part of a person's life and wellbeing. The **Aquatec Ergo Family** of shower chair commodes have been designed with user comfort, independence and dignity in mind. Next generation of the well-established Aquatec shower chairs, the Ergo family comes packed with a range of new features, designed to make personal care activities easier for both the user and carer.

Safety, independence and dignity

A newly contoured seat surface has been developed to encourage a more upright and stable seating position, promoting user independence and overall

comfort. Enhanced stability at the front of the base allows for a smaller footprint and offers user added safety and reassurance when seated in the chair. Furthermore, the optimized wheelbase makes it even easier for a carer to maneuver the chair in small spaces as well as allowing for a better fit over standard toilets for independent toileting.

Modular and flexible design

The modular design of the Ergo family means that the chairs can be adapted to suit a wide range of user needs by making some simple adjustments, or by adding some of our easy to fit, off-the-shelf accessories.

OTHER features include:

- ▶ Stable stainless steel frame
- ▶ Machine washable tension adjustable back with quick-release buckles
- ▶ Easy disassembly for transportation and storage



NEW

FEATURES

1 Innovative ergonomic seat

- Insight-led design encourages a more upright seating position for enhanced independence and comfort

2 Greater adaptability

- With an increased height adjustment range on the footrests and armrests

4 Compact and easy to handle

- Optimized wheel base makes it easy to maneuver in small spaces


3 Extremely stable

- Enhanced stability at the front for additional safety and reassurance



NEW Contoured, ergonomic seat plate

- ▶ Based on the imprints people leave behind when seated
- ▶ A 5° seat tilt lifts the knees into a natural squat position for toileting
- ▶ Transfer supports at the front to allow the user to position themselves into a comfortable position

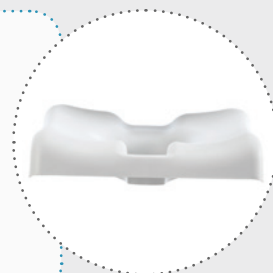

Highly Configurable

Available with a extensive range of positioning, safety, toileting and seating accessories that can be used to support a wide range of conditions and postural needs.


See full range of accessories online at www.clarkehealthcare.com



Positioning



Safety

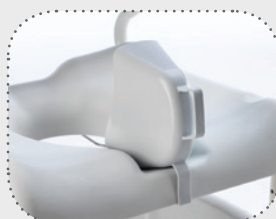

Toileting



Seating and Backrests

Tilt in space models:

Ergo VIP and Ergo Dual VIP

- ▶ Both have a -5-40° seat tilt angle with a smart weight shift mechanism that ensures optimal stability on a small footprint, whilst the Ergo Dual VIP has an additional backrest recline of 5°-35°



Introducing the **Ergo** Family

From a standard transit model, through to tilt and recline, meet the flexible range of shower chairs that offer personal care solutions for a wide range of physical conditions.



Ergo SPXL ▲

An entry level plus-size model with a maximum user weight of 400 lb.

Ergo XL ▲

An entry level plus-size attendant propelled model with a maximum user weight of 400 lb.

Ergo VIP ▲

Designed for individuals who require slightly more assistance with postural control when seated

Ergo ▲

An entry level attendant propelled model

Ergo Dual VIP ▲

Designed for individuals with complex seating and positioning needs

Ergo SP ▲

An entry level self-propelling model

Technical Data ▼

For more comprehensive information about this product, including the product's user manual, please visit www.clarkehealthcare.com.



Seat width ▼



Seat depth ▼



Total width ▼



Total depth ▼



Width between armrests ▼

	Seat width ▼	Seat depth ▼	Total width ▼	Total depth ▼	Width between armrests ▼
Ergo	18.7"	18.7"	22"	35.3"	17.75"
Ergo XL	18.7"	18.7"	25"	35.3"	20.87"
Ergo SP	18.7"	18.7"	26"	41"	17.75"
Ergo SPXL	18.7"	18.7"	26"	41"	20.87"
Ergo VIP	18.7"	18.7"	22"	39"	17.75"
Ergo Dual VIP	18.7"	18.7"	22"	39"	17.75"



Tilt angle seat ▼



Tilt angle back ▼



Weight ▼



Load capacity ▼

	Tilt angle seat ▼	Tilt angle back ▼	Weight ▼	Load capacity ▼
Ergo	-	-	37.5 lb.	330 lb.
Ergo XL	-	-	39.7 lb.	400 lb.
Ergo SP	-	-	50.7 lb.	330 lb.
Ergo SPXL	-	-	52.9 lb.	400 lb.
Ergo VIP	-5° - 40°	-	51.8 lb.	330 lb.
Ergo Dual VIP	-5° - 40°	5° - 35°	62.4 lb.	330 lb.

Frame color ▼



White

Regulatory data ▼



Clarke Healthcare Products
7830 Steubenville Pike, Oakdale, PA 15071
www.clarkehealthcare.com **888-347-4537**

