



CARROLL CS SERIES™  
**CS9™ FX600**  
 ADJUSTABLE-WIDTH BED



Yes, you can.®



CS9 FX600 bed featured with Invacare® ThinkSoft® Positioning Device, Kendall Bed Ends and Invacare® Solace® Mattress

The innovative **Carroll CS Series™ CS9™ FX600 Adjustable Width Bed** adjusts from 36" to 39" or 42" to accommodate a variety of residents, offering facilities the utmost in flexibility and providing a total solution for any long-term care facility. Matching the look and feel of the entire CS line, this FX600 is a bariatric bed that blends in with the others, bringing dignity to the resident. The optional trendelenburg option makes this bed a true solution in every wing of the facility.

Invacare beds, mattresses and rails, excluding therapeutic support surfaces and bariatric beds, when used as a system and sold new and unused after August 1, 2007, meet the voluntary FDA Guidelines for bed rail entrapment. Visit [www.invacare.com/bedsafety](http://www.invacare.com/bedsafety) for additional safety and entrapment information.



## Features

**Expandable Width Capability** – Adjusts from 36" to 39" to 42" without additional parts or tools. Tool-less adjustable headboard and 39" fixed footboard keep the dimensions within the entrapment guidelines.

**Advanced Electronics** – Synchronized motors designed to keep the bed level at all times, providing stability throughout the travel range.

**One-Touch Floor Lock** – Patented one-touch floor lock secures the bed at any height.



**Pure Vertical Rise** – Raises straight up and down throughout the travel range, with zero horizontal movement of the casters minimizing potential damage to the wall or flooring.

**Exclusive Slat Deck Design** – Designed to provide a comfortable sleep surface, allows for improved infection control and reduces wear and tear on the mattress.

**Ergonomic Hand Pendant** – Water-resistant<sup>1</sup>, rugged hand pendant designed to fit comfortably in the resident's hand with dual pendant connections.

**Customizable Options** –

- Attendant Control Panel
- Trendelenburg / Reverse Trendelenburg / Comfort Chair
- Battery Backup on board or portable
- Hi-Lo Lockout
- 4 Rail options

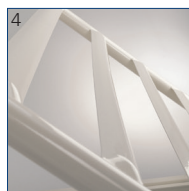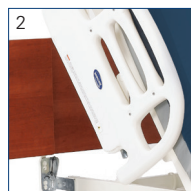

1. Expandable Deck Width
2. Expandable Headboard
3. One-Touch Floor Lock
4. Slat Deck Design
5. 6-Function Ergonomic Hand Pendant

**Invacare Corporation**  
www.invacare-cc.com

**USA**  
One Invacare Way  
Elyria, Ohio  
44035-4190  
(800) 668-2337

**Canada**  
570 Matheson Blvd. E.,  
Unit 8  
Mississauga, Ontario  
L4Z 4G4 Canada  
(800) 668-5324

©2017  
Invacare Corporation.  
All rights reserved.  
Trademarks are  
identified by the  
symbols™ and®.  
All trademarks are  
owned by or licensed  
to Invacare Corporation  
unless otherwise noted.  
Carroll CS Series and  
CS9 are trademarks  
of Carroll Healthcare.  
Specifications are  
subject to change  
without notification.  
Form No. 12-636  
Rev. 05/17  
170540



Model #	Height Range	Sleep Surface	Weight Capacity*	Regulation Compliance	Limited Warranty
IHCS9FX	9" - 30" H	76" or 80" L 36" - 39" - 42" W	600 lb.	†Classified to: UL 60601-1 IEC 60601-2-38 CAN/CSA C22.2 No. 601.1	10 Years Frame 5 Years Mechanical 3 Years Electronics

NOTE: Whenever bed rails are used as part of any hospital or homecare bed system, the potential for bed rail entrapment exists. To reduce this risk, a patient must always be assessed regarding their risk of entrapment before any product is used. Routine patient monitoring is also required once a product is in use. To learn more about the risk of bed rail entrapment please visit the US FDA's website at <http://www.fda.gov> and search for "Bed Rail Entrapment."

\* Maximum safe working load.

1. Pendant is water-resistant to IPX6 per IEC 60529.

†Testing conducted to select portions of UL/IEC Standards.

